



Instructions for Use

---


PRO-FLEX<sup>®</sup> MODULAR


		3
EN	Instructions for Use	5
DE	Gebrauchsanweisung	10
FR	Notice d'utilisation	16
ES	Instrucciones para el uso	22
IT	Istruzioni per l'uso	28
NO	Bruksanvisning	34
DA	Brugsanvisning	39
SV	Bruksanvisning	44
EL	Οδηγίες Χρήσης	49
FI	Käyttöohjeet	55
NL	Gebruiksaanwijzing	61
PT	Instruções de Utilização	67
PL	Instrukcja użytkowania	73
TR	Kullanım Talimatları	79
RU	Инструкция по использованию	84
JA	取扱説明書	90
ZH	中文说明书	95
KO	사용 설명서	100

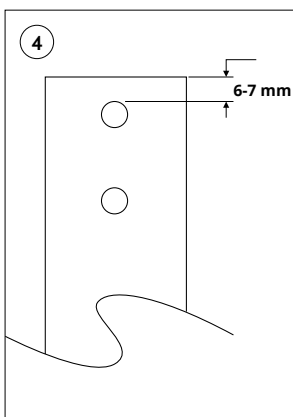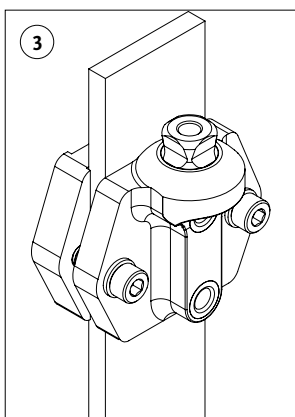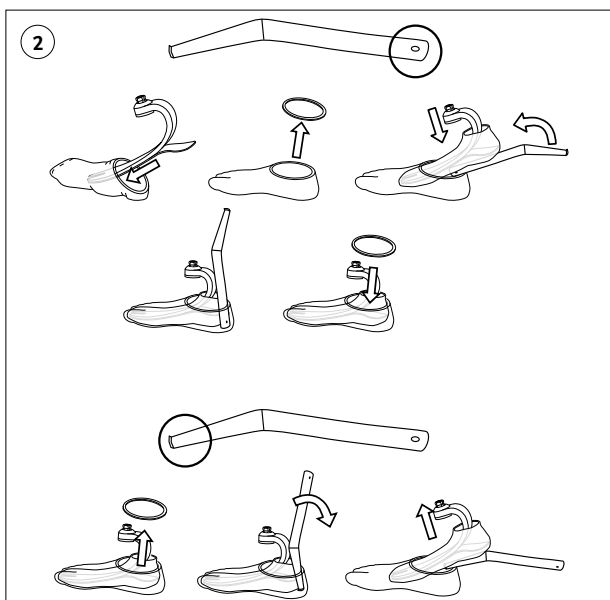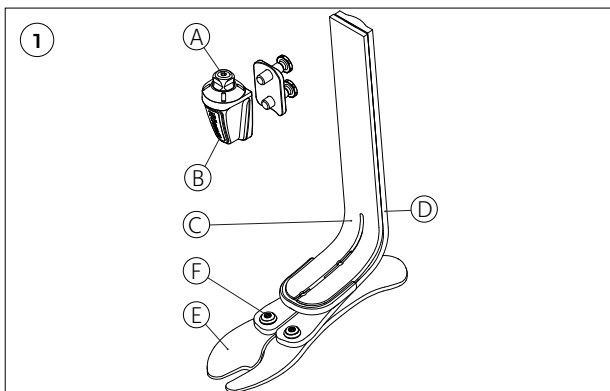5

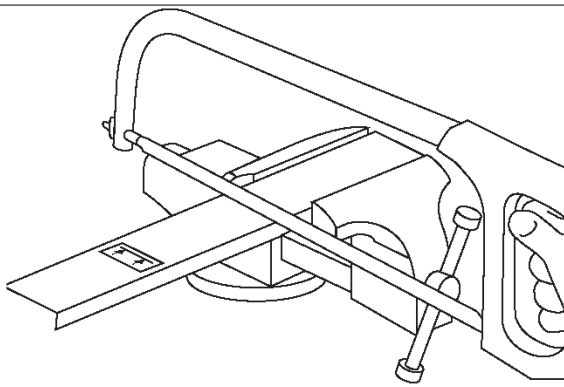6

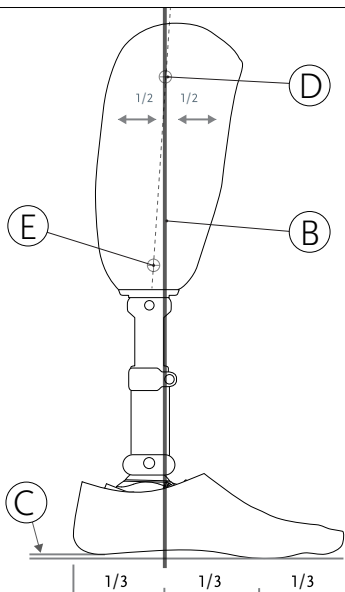Medical Device

## DESCRIPTION

The device is a prosthetic foot device with a male pyramid and three composite blades. It features a sandal toe and a full-length composite blade. The bottom blade is secured together with heel bolts.

The device consists of the following components (**Fig. 1**):

- A. Male Pyramid Adapter
- B. Pyramid Cover
- C. Top Blade
- D. Middle Blade
- E. Bottom Blade
- F. Bottom Blade Bolts

This device must be used with a Foot Cover and a Spectra Sock.

## INTENDED USE

The device is intended as a part of a prosthetic system that replaces the foot and ankle function of a missing lower limb.

Suitability of the device for the prosthesis and the patient must be evaluated by a healthcare professional.

The device must be fitted and adjusted by a healthcare professional.

### ***Indications For Use and Target Patient Population***

- Lower limb loss, amputation, or deficiency
- No known contraindications

The device is for low to high impact use, e.g., walking and occasional running.

## GENERAL SAFETY INSTRUCTIONS

**Warning:** Use of a lower limb prosthetic device carries an inherent risk of falling which may lead to injury.

The healthcare professional should inform the patient about everything in this document that is required for safe use of this device.

**Warning:** If there is a change or loss in device functionality, or if the device shows signs of damage or wear hindering its normal functions, the patient should stop using the device and contact a healthcare professional.

The device is for single patient use.

## DEVICE SELECTION

Verify that selected variant of the device is suitable for the impact level and weight limit according to the following table.

**Warning:** Do not exceed weight limit. Risk of device failure.

Incorrect category selection may also result in poor device function.

Weight kg	45-52	53-59	60-68	69-77	78-88	89-100	101-116	117-130	131-147	148-166
Weight lbs	99-115	116-130	131-150	151-170	171-194	195-220	221-256	257-287	288-324	325-365
<b>Low Impact Level</b>										
Size 22-24	1	1	2	3	4	5	6	N/A	N/A	N/A
Size 25-27	1	1	2	3	4	5	6	7	8	9
Size 28-30	3	3	3	3	4	5	6	7	8	9
<b>Moderate Impact Level</b>										
Size 22-24	1	2	3	4	5	6	N/A	N/A	N/A	N/A
Size 25-27	1	2	3	4	5	6	7	8	9	N/A
Size 28-30	3	3	3	4	5	6	7	8	9	N/A
<b>High Impact Level</b>										
Size 22-24	2	3	4	5	6	N/A	N/A	N/A	N/A	N/A
Size 25-27	2	3	4	5	6	7	8	9	N/A	N/A
Size 28-30	3	3	4	5	6	7	8	9	N/A	N/A

## ASSEMBLY INSTRUCTIONS

### **Foot Cover and Spectra Sock (Fig. 2)**

**Caution:** To avoid pinching fingers, always use a shoehorn.

1. Put the foot into the Spectra Sock.
2. Remove the foot cover attachment.
3. Use the applicable (straight) end of a shoehorn to put the foot with the Spectra Sock into the Foot Cover.
4. Move the shoehorn up to fully push the foot into the Foot Cover.
5. Put the foot cover attachment back onto the Foot Cover.
6. After alignment is complete, fix the Spectra Sock to the prosthesis to seal against dust and dirt.

**Note:** The Spectra Sock must be pulled up to prevent it from interfering with moving parts of the foot.

If required doff the Foot Cover as follows:

1. Remove the foot cover attachment.
2. Insert the applicable (edged) end of a shoehorn behind the foot.
3. Push the shoehorn down and pull the foot out of the Foot Cover.
4. Fully remove the Spectra Sock.

### **Using Modular Alignment Adapter**

#### **Alignment Adapter**

Attach the Alignment Adapter to the Foot Module (**Fig. 3**).

Connect the Foot Module and the Alignment Adapter to the prosthesis. Adjust the Alignment Adapter to the correct height by sliding it on the blade.

Clamp the Alignment Adapter to the Foot Module by tightening the screws to 22 Nm.

Disconnect the Foot Module, with the Alignment Adapter clamped onto it, from the prosthesis.

### **Drilling the Foot Module Attachment Holes**

When drilling, support the Foot Module with a wood block or a discarded carbon fiber section and clamp securely.

Use the drill guide on the Alignment Adapter to locate where to drill the holes (**Fig. 4**).

Drill 10.5 mm hole size through the drill guide.

Drill speed 1200-1800 RPM.

Feed drill slowly with light intermittent pressure, do not force.

Mark the correct length of the Foot Module so the upper edge is 6-7 mm above the top hole (**Fig. 4**).

Detach the Alignment Adapter from the Foot Module.

### **Cutting the Foot Module**

Cut the Foot Module with a sharp fine tooth hacksaw blade. Do not cut below the cut limit sticker line (**Fig. 5**).

Sand away any sharp edges.

### **Attaching the Pyramid**

Attach pyramid to the Foot Module. Place medium strength threadlocker on pyramid bolts before torquing to 45 Nm.

Connect the Foot Module, with the pyramid attached, to the prosthesis.

For secure fastening use a medium strength threadlocker and apply 15 Nm torque twice to the pyramid receiver setscrews.

### **Prosthesis**

Assemble prosthesis with applicable devices.

**Warning:** Risk of structural failure. Components from other manufacturers have not been tested and may cause excessive load on the device.

**Warning:** Ensure proper attachment by following the applicable device assembly instructions.

## **ALIGNMENT INSTRUCTIONS**

### ***Bench Alignment* (Fig. 6)**

#### **Alignment Goal**

Alignment reference line (B) should:

- pass through midpoint of socket at either the patellar tendon level or at the ischial tuberosity level (D).
- fall at the 1/3 mark on the inside of the foot cover.

**Note:** Prioritize knee alignment over foot alignment if there is a mismatch.

### **Alignment Instructions**

1. Position the foot so that the alignment reference line (B) falls at the 1/3 mark on the inside of the foot cover (with the foot cover and shoe on). Consider the external rotation of the foot.
2. Use the applicable adapters to connect either the socket or the knee to the foot and establish the correct knee center height.
3. If using a prosthetic knee: Position knee according to knee alignment instructions
4. On the lateral side of the socket, make a first mark (alignment reference point) at the midpoint of the socket at either patellar tendon level or the ischial tuberosity level (D). Make a second mark at the midpoint of the socket distally (E). Draw a line through both marks.
5. Position the socket so the alignment reference line (B) goes through the alignment reference point (D).
6. Set the correct socket angles for flexion/extension and abduction/adduction.
7. If using a prosthetic knee: Use the applicable adapters to connect the socket to the prosthesis.

### ***Static Alignment***

- Make sure the patient stands with equal weight on both legs.
- Check for correct prosthesis length.

- Check internal / external rotation.
- Check for correct load on toe and heel.

### ***Dynamic Alignment***

Ensure the patient is familiar with the functioning of the device.

The heel to toe action can be influenced by:

- Heel Stiffness
- Anterior-Posterior positioning of device
- Dorsi-Plantarflexion
- Shoe Characteristics

Consider the following actions if needed:

#### **Symptoms**

- Device comes to flat position too early (patient feels they are sinking into a hole)
- Climbing over the toe requires extra energy
- Toe feels too stiff
- Knee hyperextends

#### **Action**

- Check Wedge instructions (if applicable)
- Shift socket anterior (or device posterior)
- Consider dorsiflexion
- Check heel of the shoe and shoe performance

#### **Symptoms**

- Rapid heel to toe movement
- Poor control over prosthesis at initial contact
- Minimal energy return feeling
- Too little push off from the toe
- Knee becomes unstable

#### **Action**

- Check Wedge instructions (if applicable)
- Shift socket posterior (or device anterior)
- Consider plantarflexion
- Check heel of the shoe and shoe performance

## **USAGE**

### ***Cleaning and care***

Clean with a damp cloth and a mild soap. Dry with a cloth after cleaning.

### ***Environmental Conditions***

The device is Waterproof.

A Waterproof device can be used in a wet or humid environment and submerged in up to 3-meter-deep water for a maximum of 1 hour.

It can tolerate contact with: Salt water, chlorinated water, perspiration, urine, and mild soaps.

It can also tolerate occasional exposure to sand, dust, and dirt.

Continuous exposure is not allowed.

Dry with a cloth after contact with fresh water or humidity.

Clean with fresh water and dry with a cloth after exposure to other liquids and chemicals or in case of accidental exposure to sand, dust, or dirt.

## **MAINTENANCE**

The device and the overall prosthesis should be examined by a healthcare professional. Interval should be determined based on patient activity.

### Noise from Foot

Noise may occur if sand or debris is present in device. In that case, the healthcare professional should doff the foot, clean it with the help of compressed air and replace the Spectra sock if it is damaged.

### REPORT OF SERIOUS INCIDENT

Any serious incident in relation to the device must be reported to the manufacturer and relevant authorities.

### DISPOSAL

The device and packaging must be disposed of in accordance with respective local or national environmental regulations.

### LIABILITY

Össur prosthetic devices are designed and verified to be safe and compatible in combination with each other and custom-made prosthetic sockets with Össur adapters, and when used in accordance with their intended use.

Össur does not assume liability for the following:

- Device not maintained as instructed by the instructions for use.
- Device assembled with components from other manufacturers.
- Device used outside of recommended use condition, application, or environment.

### Compliance

This device has been tested according to ISO 10328 standard to two million load cycles.

Depending on patient activity this may correspond to 2-3 years of use.

ISO 10328 - label		
Category	Weight (Kg) Moderate Impact	Label Text
1	52	ISO 10328-P3*-52 kg
2	59	ISO 10328-P3*-59 kg
3	68	ISO 10328-P3-68 kg
4	77	ISO 10328-P3-77 kg
5	88	ISO 10328-P4-88 kg
6	100	ISO 10328-P5-100 kg
7	116	ISO 10328-P6-116 kg
8	130	ISO 10328-P7-130 kg
9	147	ISO 10328-P8-147 kg
10	166	ISO 10328-P8-166 kg
11	185	ISO 10328-P8-185 kg
12	203	ISO 10328-P8-203 kg

\*Actual test loads reflect maximum body mass

ISO 10328 - "P" - "m"kg \*)

\*) Body mass limit not to be exceeded!

For specific conditions and limitations of use see manufacturer's written instructions on intended use!

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