



Instructions for Use

---

PRO-FLEX<sup>®</sup> TERRA

		3
EN	Instructions for Use	4
DE	Gebrauchsanweisung	9
FR	Notice d'utilisation	14
ES	Instrucciones para el uso	19
IT	Istruzioni per l'uso	24
NO	Bruksanvisning	29
DA	Brugsanvisning	34
SV	Bruksanvisning	39
EL	Οδηγίες Χρήσης	44
FI	Käyttöohjeet	49
NL	Gebruiksaanwijzing	54
PT	Instruções de Utilização	59
PL	Instrukcja użytkowania	64
TR	Kullanım Talimatları	69
RU	Инструкция по использованию	74
JA	取扱説明書	79
ZH	中文说明书	84
KO	사용 설명서	89
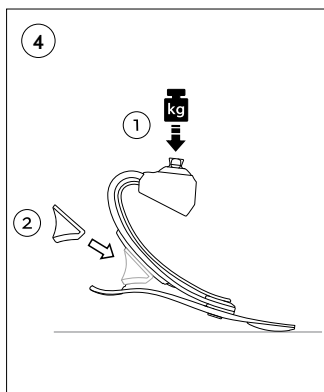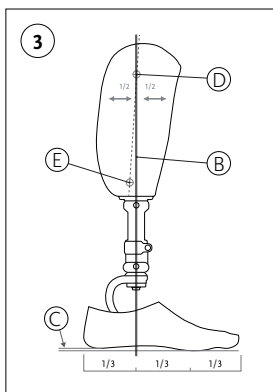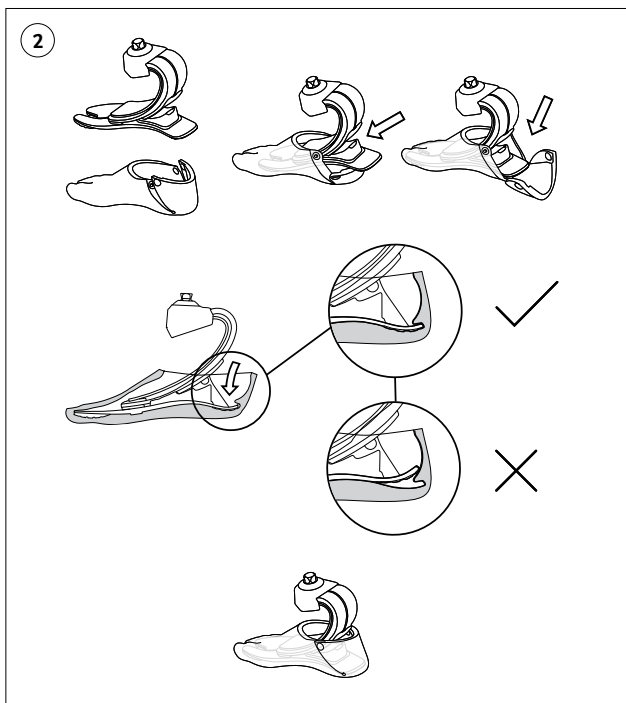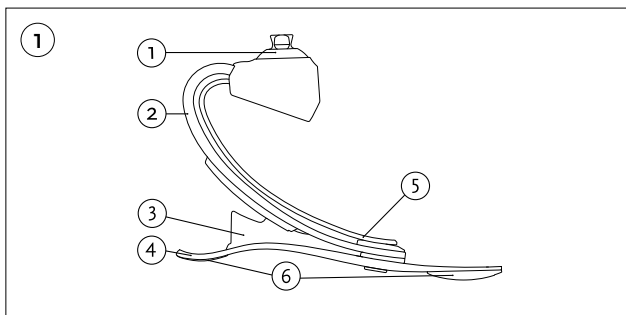
### DESCRIPTION

The device is a multi-blade prosthetic foot featuring a pre-compressed top blade.

The device consists of the following components (**Fig. 1**):

- A. Pyramid
- B. Middle blade
- C. Heel foam
- D. Bottom blade
- E. Top blade
- F. Grip pads

The device must be used with a foot cover.

### INTENDED USE

The device is intended as a part of a prosthetic system that replaces the foot and ankle function of a missing lower limb.

Suitability of the device for the prosthesis and the patient must be evaluated by a healthcare professional.

The device must be fitted and adjusted by a healthcare professional.

### *Indications For Use and Target Patient Population*

- Lower limb loss, amputation, or deficiency
- No known contraindications

The device is for low to high impact use, e.g., walking and occasional running.

### GENERAL SAFETY INSTRUCTIONS

**Warning:** Use of a lower limb prosthetic device carries an inherent risk of falling which may lead to injury.

The healthcare professional should inform the patient about everything in this document that is required for safe use of this device.

**Warning:** If there is a change or loss in device functionality, or if the device shows signs of damage or wear hindering its normal functions, the patient should stop using the device and contact a healthcare professional.

The device is for single patient use.

### DEVICE SELECTION

Verify that selected variant of the device is suitable for the impact level and weight limit according to the following table.

**Warning:** Do not exceed weight limit. Risk of device failure.

Incorrect category selection may also result in poor device function.

Weight kg	45-52	53-59	60-68	69-77	78-88	89-100	101-116	117-130	131-147	148-166
Weight lbs	99-115	116-130	131-150	151-170	171-194	195-220	221-256	257-287	288-324	325-365
<b>Low Impact Level</b>										
Size 22-24	1	1	2	3	4	5	6	N/A	N/A	N/A
Size 25-27	1	1	2	3	4	5	6	7	8	9
Size 28-30	3	3	3	3	4	5	6	7	8	9
<b>Moderate Impact Level</b>										
Size 22-24	1	2	3	4	5	6	N/A	N/A	N/A	N/A
Size 25-27	1	2	3	4	5	6	7	8	9	N/A
Size 28-30	3	3	3	4	5	6	7	8	9	N/A
<b>High Impact Level</b>										
Size 22-24	2	3	4	5	6	N/A	N/A	N/A	N/A	N/A
Size 25-27	2	3	4	5	6	7	8	9	N/A	N/A
Size 28-30	3	3	4	5	6	7	8	9	N/A	N/A

## ASSEMBLY INSTRUCTIONS

### Foot Cover (Fig. 2)

**Warning:** Using the device without a foot cover may result in a fall.

To don the foot cover:

1. Open the snaps and pull the heel section down.
2. Insert the foot into the foot cover.
3. Make sure the bottom blade end is inside the heel groove.
4. Close both snaps.

To doff the foot cover:

1. Open the snaps and pull the heel section down.
2. Remove the foot.

### Prosthesis

Assemble prosthesis with applicable devices.

**Warning:** Risk of structural failure. Components from other manufacturers have not been tested and may cause excessive load on the device.

**Warning:** Ensure proper attachment by following the applicable device assembly instructions.

Remove the Pyramid Protector during final fitting.

**Caution:** The Pyramid Protector may have sharp edges.

## ALIGNMENT INSTRUCTIONS

### Bench Alignment (Fig. 3)

#### Alignment Goal

Alignment reference line (B) should:

- pass through midpoint of socket at either the patellar tendon level or at the ischial tuberosity level (D).
- fall at the 1/3 mark on the inside of the foot cover.

**Note:** Prioritize knee alignment over foot alignment if there is a mismatch.

#### Alignment Instructions

1. Position the foot so that the alignment reference line (B) falls at the 1/3 mark on the inside of the foot cover (with the foot cover and shoe on). Consider the external rotation of the foot.

2. Use the applicable adapters to connect either the socket or the knee to the foot and establish the correct knee center height.
3. If using a prosthetic knee: Position knee according to knee alignment instructions
4. On the lateral side of the socket, make a first mark (alignment reference point) at the midpoint of the socket at either patellar tendon level or the ischial tuberosity level (D). Make a second mark at the midpoint of the socket distally (E). Draw a line through both marks.
5. Position the socket so the alignment reference line (B) goes through the alignment reference point (D).
6. Set the correct socket angles for flexion/extension and abduction/adduction.
7. If using a prosthetic knee: Use the applicable adapters to connect the socket to the prosthesis.

### **Static Alignment**

- Make sure the patient stands with equal weight on both legs.
- Check for correct prosthesis length.
- Check internal / external rotation.
- Check for correct load on toe and heel.

### **Dynamic Alignment**

Ensure the patient is familiar with the functioning of the device.

The heel to toe action can be influenced by:

- Heel Stiffness
- Anterior-Posterior positioning of device
- Dorsi-Plantarflexion
- Shoe Characteristics

Consider the following actions if needed:

### **Symptoms**

- Device comes to flat position too early (patient feels he/she is sinking into a hole)
- Climbing over the toe requires extra energy
- Toe feels too stiff
- Knee hyperextends

### **Action**

- Check Wedge instructions (if applicable)
- Shift socket anterior (or device posterior)
- Consider dorsiflexion
- Check heel of the shoe and shoe performance

### **Symptoms**

- Rapid heel to toe movement
- Poor control over prosthesis at initial contact
- Minimal energy return feeling
- Too little push off from the toe
- Knee becomes unstable

### **Action**

- Check Wedge instructions (if applicable)
- Shift socket posterior (or device anterior)
- Consider plantarflexion
- Check heel of the shoe and shoe performance

### ***Heel Wedge***

Heel wedges are used to increase heel stiffness. They can be used for short or long time.

To install heel wedge, load the toe of the foot and push wedge in place.

Make sure that the wedge tab fits into the groove of the Heel foam (**Fig. 4**).

To remove heel wedge, load the toe of the foot and pull wedge out.

## **USAGE**

### ***Cleaning and care***

Remove foot cover to clean. Rinse with water and remove all sand, dust, and dirt from foot and foot cover. Dry with a cloth after cleaning.

### ***Environmental Conditions***

The device is Waterproof.

A Waterproof device can be used in a wet or humid environment and submerged in up to 3-meter-deep water for a maximum of 1 hour.

It can tolerate contact with: Salt water, chlorinated water, perspiration, urine, and mild soaps.

It can also tolerate occasional exposure to sand, dust, and dirt.

Continuous exposure is not allowed.

Dry with a cloth after contact with fresh water or humidity.

Clean with fresh water after exposure to other liquids, chemicals, sand, dust, or dirt and dry with a cloth.

## **MAINTENANCE**

The device and the overall prosthesis should be examined by a healthcare professional. Interval should be determined based on patient activity.

### ***Noise from Foot***

Noise may occur if sand or debris is present in device. In that case, clean device and foot cover. If noise persists, a Spectra sock can be used to eliminate noise between device and foot cover.

## **REPORT OF SERIOUS INCIDENT**

Any serious incident in relation to the device must be reported to the manufacturer and relevant authorities.

## **DISPOSAL**

The device and packaging must be disposed of in accordance with respective local or national environmental regulations.

## **LIABILITY**

Össur prosthetic devices are designed and verified to be safe and compatible in combination with each other and custom-made prosthetic sockets with Össur adapters, and when used in accordance with their intended use.

Össur does not assume liability for the following:

- Device not maintained as instructed by the instructions for use.
- Device assembled with components from other manufacturers.
- Device used outside of recommended use condition, application, or environment.

## Compliance

This device has been tested according to ISO 10328 standard to two million load cycles.

Depending on patient activity this may correspond to 2-3 years of use.

ISO 10328 - Label		
Category	Weight (Kg) Moderate Impact	Label Text
1	52	ISO 10328-P3*-52 kg
2	59	ISO 10328-P3*-59 kg
3	68	ISO 10328-P3-68 kg
4	77	ISO 10328-P3-77 kg
5	88	ISO 10328-P4-88 kg
6	100	ISO 10328-P5-100 kg
7	116	ISO 10328-P6-116 kg
8	130	ISO 10328-P7-130 kg
9	147	ISO 10328-P8-147 kg

\*Actual test loads reflect maximum body mass

ISO 10328 - "P" - "m"kg \*)

\*) Body mass limit not to be exceeded!

For specific conditions and limitations of use see manufacturer's written instructions on intended use!

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