










Instructions for Use

RE-FLEX ROTATE™
RE-FLEX SHOCK™

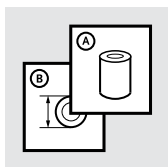		3
EN	Instructions for Use	5
DE	Gebrauchsanweisung	12
FR	Notice d'utilisation	20
ES	Instrucciones para el uso	28
IT	Istruzioni per l'uso	35
NO	Bruksanvisning	42
DA	Brugsanvisning	49
SV	Bruksanvisning	56
EL	Οδηγίες Χρήσης	63
FI	Käyttöohjeet	71
NL	Gebruiksaanwijzing	78
PT	Instruções de Utilização	86
PL	Instrukcja użytkowania	93
TR	Kullanım Talimatları	100
RU	Инструкция по использованию	107
JA	取扱説明書	115
ZH	中文说明书	122
KO	사용 설명서	129

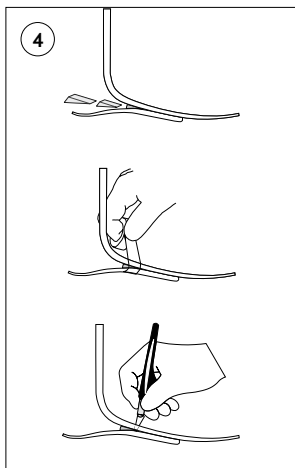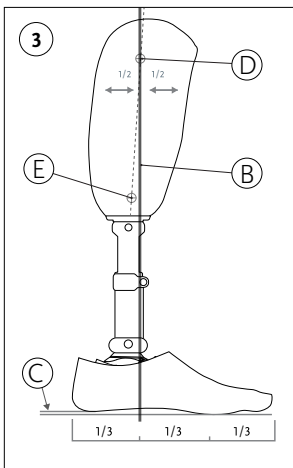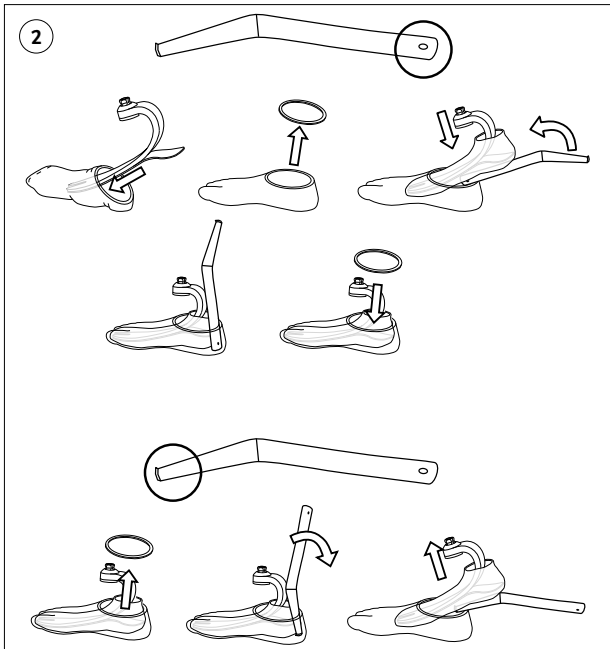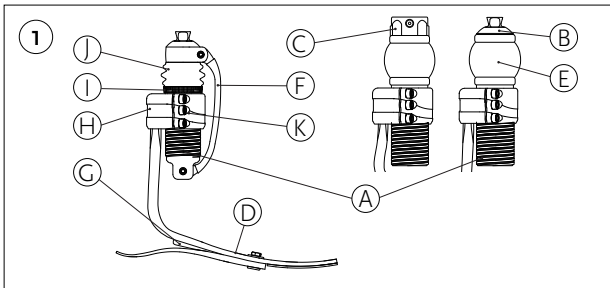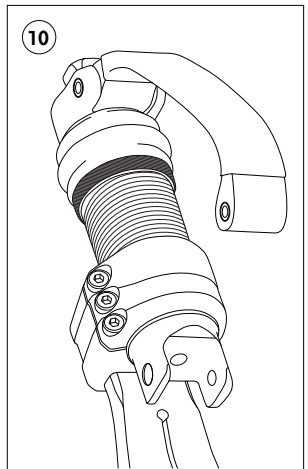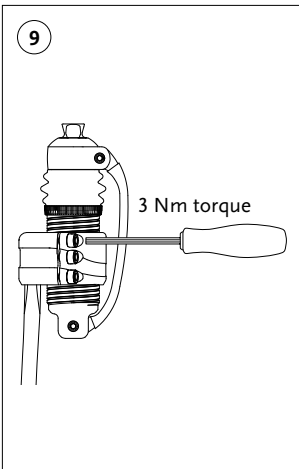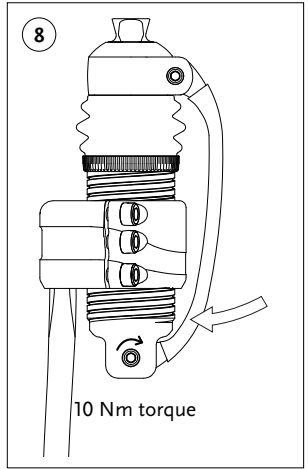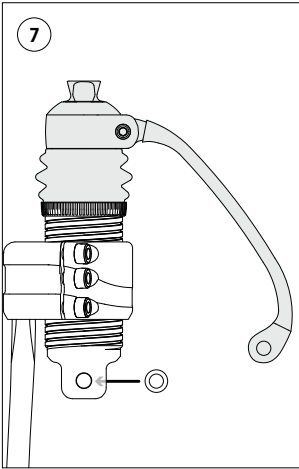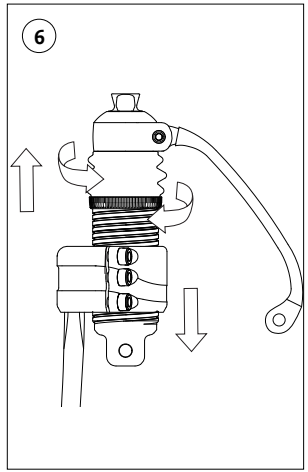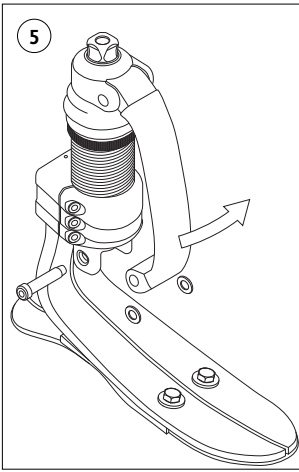DESCRIPTION

Re-Flex Rotate is a mechanical foot with an integrate torsion and vertical shock absorbing adapter, available with male or female pyramid options. Re-Flex Shock is a mechanical foot with an integrated vertical shock-absorbing adapter and a male pyramid.

The device consists of the following components (**Fig. 1**):

- A. Shock Module
- B. Male Pyramid
- C. Female Pyramid
- D. Foot Module
- E. Torsion Cell
- F. Front Spring
- G. Heel Wedge
- H. Attachment Clamp
- I. Collar
- J. Bellow

This device must be used with a Foot Cover and a Spectra Sock.

INTENDED USE

The device is intended as a part of a prosthetic system that replaces the foot and ankle function of a missing lower limb.

Suitability of the device for the prosthesis and the patient must be evaluated by a healthcare professional.

The device must be fitted and adjusted by a healthcare professional.

Indications For Use and Target Patient Population

- Lower limb loss, amputation, or deficiency
- No known contraindications

The device is for low to high impact use, e.g., walking and occasional running.

GENERAL SAFETY INSTRUCTIONS

Warning: Use of a lower limb prosthetic device carries an inherent risk of falling which may lead to injury.

The healthcare professional should inform the patient about everything in this document that is required for safe use of this device.

Warning: If there is a change or loss in device functionality, or if the device shows signs of damage or wear hindering its normal functions, the patient should stop using the device and contact a healthcare professional.

The device is for single patient use.

DEVICE SELECTION

Verify that selected variant of the device is suitable for the impact level and weight limit according to the following table.

Warning: Do not exceed weight limit. Risk of device failure.

Incorrect category selection may also result in poor device function.

Re-Flex Rotate

Weight kg	45-52	53-59	60-68	69-77	78-88	89-100	101-116	117-130	131-147
Weight lbs	99-115	116-130	131-150	151-170	171-194	195-220	221-256	257-287	288-324
Low Impact Level									
Size 22-24	1	1	2	3	4	5	6	N/A	N/A
Size 25-27	1	1	2	3	4	5	6	7	8
Size 28-30	3	3	3	3	4	5	6	7	8
Moderate Impact Level									
Size 22-24	1	2	3	4	5	6	N/A	N/A	N/A
Size 25-27	1	2	3	4	5	6	7	8	N/A
Size 28-30	3	3	3	4	5	6	7	8	N/A
High Impact Level									
Size 22-24	2	3	4	5	N/A	N/A	N/A	N/A	N/A
Size 25-27	2	3	4	5	6	7	8	N/A	N/A
Size 28-30	3	3	4	5	6	7	8	N/A	N/A

Re-Flex Shock

Weight kg	45-52	53-59	60-68	69-77	78-88	89-100	101-116	117-130	131-147	148 - 166
Weight lbs	99-115	116-130	131-150	151-170	171-194	195-220	221-256	257-287	288-324	325-366
Low Impact Level										
Size 22-24	1	1	2	3	4	5	6	N/A	N/A	N/A
Size 25-27	1	1	2	3	4	5	6	7	8	9
Size 28-30	3	3	3	3	4	5	6	7	8	9
Moderate Impact Level										
Size 22-24	1	2	3	4	5	6	N/A	N/A	N/A	N/A
Size 25-27	1	2	3	4	5	6	7	8	9	N/A
Size 28-30	3	3	3	4	5	6	7	8	9	N/A
High Impact Level										
Size 22-24	2	3	4	5	6	N/A	N/A	N/A	N/A	N/A
Size 25-27	2	3	4	5	6	7	8	9	N/A	N/A
Size 28-30	3	3	4	5	6	7	8	9	N/A	N/A

ASSEMBLY INSTRUCTIONS

Foot Cover and Spectra Sock (Fig. 2)

Caution: To avoid pinching fingers, always use a shoehorn.

1. Put the foot into the Spectra Sock.
2. Remove the foot cover attachment.
3. Use the applicable (straight) end of a shoehorn to put the foot with the Spectra Sock into the Foot Cover.
4. Move the shoehorn up to fully push the foot into the Foot Cover.
5. Put the foot cover attachment back onto the Foot Cover.
6. After alignment is complete, fix the Spectra Sock to the prosthesis to seal against dust and dirt.

Note: The Spectra Sock must be pulled up to prevent it from interfering with moving parts of the foot.

If required doff the Foot Cover as follows:

1. Remove the foot cover attachment.
2. Insert the applicable (edged) end of a shoehorn behind the foot.
3. Push the shoehorn down and pull the foot out of the Foot Cover.
4. Fully remove the Spectra Sock.

Prosthesis

Assemble prosthesis with applicable devices.

Warning: Risk of structural failure. Components from other manufacturers have not been tested and may cause excessive load on the device.

Warning: Ensure proper attachment by following the applicable device assembly instructions.

Remove the protective film on the pyramid after fitting.

ALIGNMENT INSTRUCTIONS

Bench Alignment (Fig. 3)

Alignment Goal

Alignment reference line (B) should:

- pass through midpoint of socket at either the patellar tendon level or at the ischial tuberosity level (D).
- fall at the 1/3 mark on the inside of the foot cover.

Note: Prioritize knee alignment over foot alignment if there is a mismatch.

Alignment Instructions

1. Position the foot so that the alignment reference line (B) falls at the 1/3 mark on the inside of the foot cover (with the foot cover and shoe on). Consider the external rotation of the foot.
2. Use the applicable adapters to connect either the socket or the knee to the foot and establish the correct knee center height.
3. If using a prosthetic knee: Position knee according to knee alignment instructions
4. On the lateral side of the socket, make a first mark (alignment reference point) at the midpoint of the socket at either patellar tendon level or the ischial tuberosity level (D). Make a second mark at the midpoint of the socket distally (E). Draw a line through both marks.
5. Position the socket so the alignment reference line (B) goes through the alignment reference point (D).
6. Set the correct socket angles for flexion/extension and abduction/adduction.
7. If using a prosthetic knee: Use the applicable adapters to connect the socket to the prosthesis.

Static Alignment

- Make sure the patient stands with equal weight on both legs.
- Check for correct prosthesis length.
- Check internal / external rotation.
- Check for correct load on toe and heel.

Dynamic Alignment

Ensure the patient is familiar with the functioning of the device.

The heel to toe action can be influenced by:

- Heel Stiffness
- Anterior-Posterior positioning of device
- Dorsi-Plantarflexion
- Shoe Characteristics

Consider the following actions if needed:

Symptoms

- Device comes to flat position too early (patient feels he/she is sinking into a hole)

- Climbing over the toe requires extra energy
- Toe feels too stiff
- Knee hyperextends

Action

- Check Wedge instructions (if applicable)
- Shift socket anterior (or device posterior)
- Consider dorsiflexion
- Check heel of the shoe and shoe performance

Symptoms

- Rapid heel to toe movement
- Poor control over prosthesis at initial contact
- Minimal energy return feeling
- Too little push off from the toe
- Knee becomes unstable

Action

- Check Wedge instructions (if applicable)
- Shift socket posterior (or device anterior)
- Consider plantarflexion
- Check heel of the shoe and shoe performance

Heel Wedges (Fig. 4)

The small, medium and large Wedges are used to change heel stiffness. The Wedges can be trimmed using sharp scissors to customize stiffness. For temporary Wedge placement, use tape to secure the wedge in position. For permanent Wedge placement

- Roughen the upper and lower surface of the heel with abrasive paper.
- Apply instant adhesive on the lower side of the Wedge only.
- Locate in the foot/heel junction and position before adhesive sets.
- If the foot has a split toe, remove a thin slice in the middle by cutting with a sharp knife through the split in the carbon blades.

For removal, the adhesive may be softened by soaking in acetone or cyanoacrylate adhesive remover.

Re-Flex Rotate

1. Refer to **Bench Alignment**
2. **Adjusting the Shock Module**

To adjust the Shock Module, loosen the three hex screws on the Attachment Clamp (**Fig. 6**) and rotate the foot until the correct height and/or rotation is achieved. One full rotation results in 8 mm change in height. For finer adjustments, rotate in 90° turns.

Warning: Ensure there is full engagement of the clamp and Shock Module.

3. **Securing the Shock Module**

To secure the Shock Module tighten the 3 hex screws (**Fig. 9**). Start with the top hex screw and work down by applying 3 Nm. Perform this sequence twice since the upper hex screws will loosen in the first run.

Warning: Overtightening the hex screws can strip threads and cause Shock Module to bind.

4. **Applying a Permanent Shock Module Setting**

To permanently fix the Shock Module, remove the hex screws one by one; apply a medium strength threadlocker and torque of 3 Nm. When finished, re-apply torque of 3 Nm on all hex screws. Failing to apply threadlocker to the hex screws may result in the shock coming loose.

Re-Flex Shock

1. Refer to **Bench Alignment**
2. **Adjusting the Shock Module**

To adjust the Shock Module, loosen lower pin on the spring (**Fig. 5**) and the three hex screws (**Fig. 9**) on the Attachment Clamp. Rotate the shock housing until correct height is achieved (**Fig. 6**). One full rotation results in 8 mm change in height, height can be set in half-turn steps of 4 mm. The Front Spring (shown in grey color in **Fig. 7**) can be rotated 180° in relation to the threaded tube (shown in white) so the spring can always be positioned at the front when adjusting the height.

3. **Securing the Shock Module**

To secure the Shock Module fasten the lower pin on the spring (**Fig. 8**), place one disc spring on either side (**Fig. 7**), applying 10 Nm and tighten the 3 hex screws (**Fig. 9**). Start with the top hex screw and work down by applying 3 Nm. Perform this sequence twice, since the upper hex screw will loosen in the first run. Place the Bellow in the groove on the threaded tube and then tighten the Collar over the Bellow. (**Fig. 10**)

Note: Failing to add a disc spring on each side of the front spring will result in clicking noise of the foot during unloading.

4. **Applying a Permanent Shock Module Setting**

To permanently fix the Shock Module, remove the hex screws one by one; apply a medium strength threadlocker. When finished, re-apply torque of 3 Nm on all hex screws. Perform this sequence twice, since the upper hex screw will loosen in the first run. Failing to apply threadlocker to the hex screws may result in the shock coming loose. Also remove the pin and apply threadlocker on threads and tighten with 10 Nm torque.

Note: Make sure the Collar covers the edge of the Bellow all around to avoid dust and debris to come into the Shock Module.

Warning: Over tightening the hex screws can strip threads and cause the Shock Module to bind.

Re-Flex Rotate and Re-Flex Shock

Ensure there is full engagement of the clamp and Shock Module. Leave a small gap between the Collar and the clamp.

USAGE

Cleaning and care

Clean with a damp cloth and a mild soap. Dry with a cloth after cleaning. The Shock Module is lifetime lubricated and should NOT be lubricated.

Clicking Noise

If clicking noise occurs:

Check that the adapter screws are torqued and secured with the recommended values and no debris is in the interconnecting zone of the pyramid dome.

Check or remove dust between the heel and the foot blade. Clean and / or replace the Foot Cover as well as the Spectra Sock.

For Re-Flex Shock specific, check that the disc springs are inserted into the gap between the lower section of the outer tube and the Front Spring (**Fig. 7**).

Environmental Conditions

The device is Weatherproof.

A Weatherproof device can be used in a wet or humid environment and can tolerate being splashed by fresh water (e. g., rain), no submersion is allowed.

No contact with salt water or chlorinated water is allowed.

Dry with a cloth after contact with fresh water or humidity.

Clean with fresh water in case of accidental exposure to other liquids, chemicals, sand, dust, or dirt and dry with a cloth.

MAINTENANCE

The device and the overall prosthesis should be examined by a healthcare professional. Interval should be determined based on patient activity.

Noise from Foot

Noise may occur if sand or debris is present in device. In that case, the healthcare professional should doff the foot, clean it with the help of compressed air and replace the Spectra sock if it is damaged.

REPORT OF SERIOUS INCIDENT

Any serious incident in relation to the device must be reported to the manufacturer and relevant authorities.

DISPOSAL

The device and packaging must be disposed of in accordance with respective local or national environmental regulations.

LIABILITY

Össur prosthetic devices are designed and verified to be safe and compatible in combination with each other and custom-made prosthetic sockets with Össur adapters, and when used in accordance with their intended use.

Össur does not assume liability for the following:

- Device not maintained as instructed by the instructions for use.
- Device assembled with components from other manufacturers.
- Device used outside of recommended use condition, application, or environment.

Compliance

This device has been tested according to ISO 10328 standard to two million load cycles.

Depending on patient activity this may correspond to 2-3 years of use.

Re-Flex Rotate

ISO 10328 - label		
Category	Weight (Kg) Moderate Impact	Label Text
1	52	ISO 10328-P3*-52 kg
2	59	ISO 10328-P3*-59 kg
3	68	ISO 10328-P3-68 kg
4	77	ISO 10328-P3-77 kg
5	88	ISO 10328-P4-88 kg
6	100	ISO 10328-P5-100 kg
7	116	ISO 10328-P6-116 kg
8	130	ISO 10328-P7-130 kg

*Actual test loads reflect maximum body mass

Re-Flex Shock

ISO 10328 - label		
Category	Weight (Kg) Moderate Impact	Label Text
1	52	ISO 10328-P3*-52 kg
2	59	ISO 10328-P3*-59 kg
3	68	ISO 10328-P3-68 kg
4	77	ISO 10328-P3-77 kg
5	88	ISO 10328-P4-88 kg
6	100	ISO 10328-P5-100 kg
7	116	ISO 10328-P6-116 kg
8	130	ISO 10328-P7-130 kg
9	147	ISO 10328-P8-147 kg

*Actual test loads reflect maximum body mass

ISO 10328 - "P" - "m"kg *)

*) Body mass limit not to be exceeded!

For specific conditions and limitations of use
see manufacturer's written instructions on
intended use!

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