
## Instructions for Use

---


VARI-FLEX®

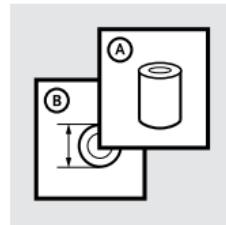
3

EN   Instructions for Use	4
DE   Gebrauchsanweisung	9
FR   Notice d'utilisation	14
ES   Instrucciones para el uso	19
IT   Istruzioni per l'uso	24
NO   Bruksanvisning	29
DA   Brugsanvisning	34
SV   Bruksanvisning	39
EL   Οδηγίες Χρήσης	44
FI   Käyttöohjeet	49
NL   Gebruiksaanwijzing	54
PT   Instruções de Utilização	59
PL   Instrukcja użytkowania	64
TR   Kullanım Talimatları	69
RU   Инструкция по использованию	74
JA   取扱説明書	79
ZH   中文说明书	84
KO   사용 설명서	89


1

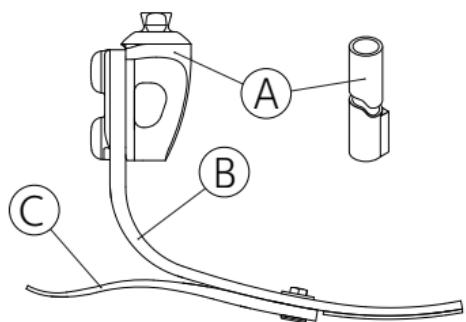
2

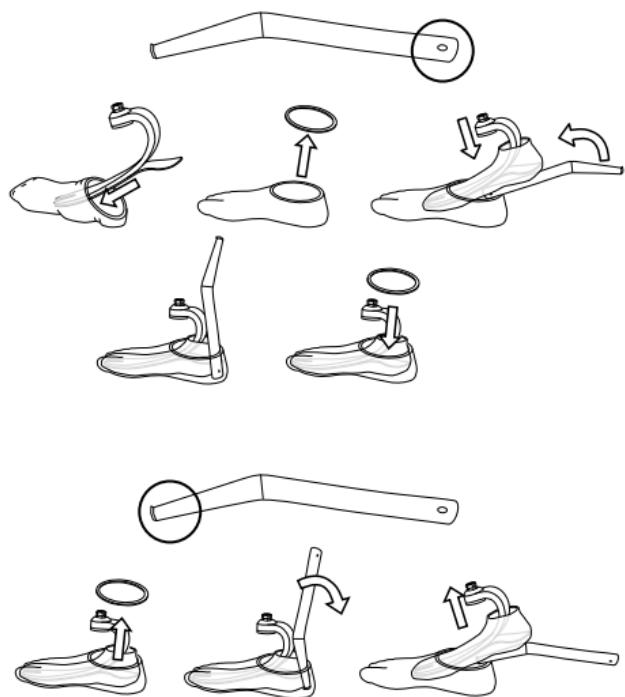3

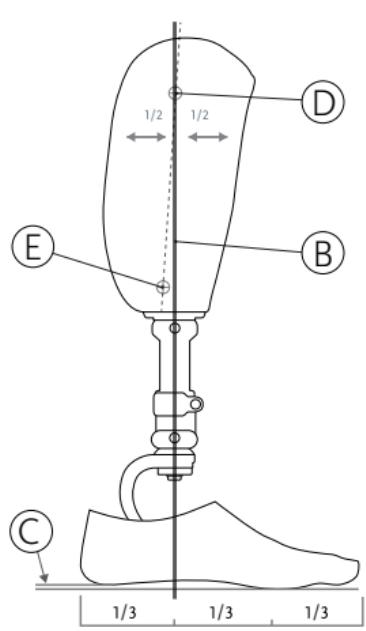4
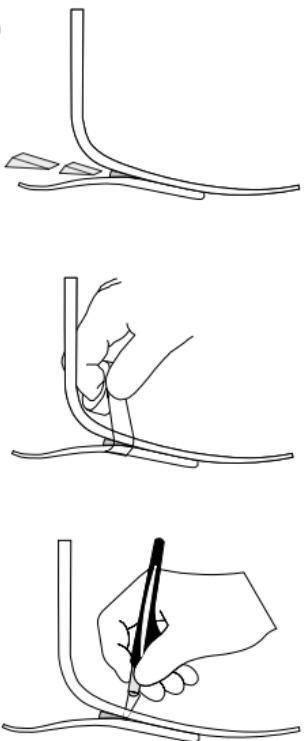
# ENGLISH

---

Medical Device

## DESCRIPTION

The device is a prosthetic foot. It is composed of a base foot, with a carbon-fiber heel and keel, and male pyramid.

The device consists of the following components (**Fig. 1**):

A. Male Pyramid or Carbon Tube

B. Foot Blade

C. Heel

This device must be used with a Foot Cover and a Spectra Sock.

## INTENDED USE

The device is intended as a part of a prosthetic system that replaces the foot and ankle function of a missing lower limb.

Suitability of the device for the prosthesis and the patient must be evaluated by a healthcare professional.

The device must be fitted and adjusted by a healthcare professional.

### ***Indications For Use and Target Patient Population***

- Lower limb loss, amputation, or deficiency
- No known contraindications

The device is for low to high impact use, e.g., walking and occasional running.

## GENERAL SAFETY INSTRUCTIONS

**Warning:** Use of a lower limb prosthetic device carries an inherent risk of falling which may lead to injury.

The healthcare professional should inform the patient about everything in this document that is required for safe use of this device.

**Warning:** If there is a change or loss in device functionality, or if the device shows signs of damage or wear hindering its normal functions, the patient should stop using the device and contact a healthcare professional.

The device is for single patient use.

## DEVICE SELECTION

Verify that selected variant of the device is suitable for the impact level and weight limit according to the following table.

**Warning:** Do not exceed weight limit. Risk of device failure.

Incorrect category selection may also result in poor device function.

Weight kg	45-52	53-59	60-68	69-77	78-88	89-100	101-116	117-130	131-147	148-166
Weight lbs	99-115	116-130	131-150	151-170	171-194	195-220	221-256	257-287	288-324	325-365
Low Impact Level	1	1	2	3	4	5	6	7	8	9
Moderate Impact Level	1	2	3	4	5	6	7	8	9	N/A
High Impact Level	2	3	4	5	6	7	8	9	N/A	N/A

## ASSEMBLY INSTRUCTIONS

### ***Foot Cover and Spectra Sock (Fig. 2)***

**Caution:** To avoid pinching fingers, always use a shoehorn.

1. Put the foot into the Spectra Sock.
2. Remove the foot cover attachment.
3. Use the applicable (straight) end of a shoehorn to put the foot with the Spectra Sock into the Foot Cover.
4. Move the shoehorn up to fully push the foot into the Foot Cover.
5. Put the foot cover attachment back onto the Foot Cover.
6. After alignment is complete, fix the Spectra Sock to the prosthesis to seal against dust and dirt.

**Note:** The Spectra Sock must be pulled up to prevent it from interfering with moving parts of the foot.

If required doff the Foot Cover as follows:

1. Remove the foot cover attachment.
2. Insert the applicable (edged) end of a shoehorn behind the foot.
3. Push the shoehorn down and pull the foot out of the Foot Cover.
4. Fully remove the Spectra Sock.

### ***Prosthesis***

Assemble prosthesis with applicable devices.

**Warning:** Risk of structural failure. Components from other manufacturers have not been tested and may cause excessive load on the device.

**Warning:** Ensure proper attachment by following the applicable device assembly instructions.

Remove the protective film on the pyramid after fitting.

### ***Bolt Assembly***

To attach a Pyramid or a Tube, thread each bolt through washer and Foot Module into Pyramid/Tube. Apply Loctite 410 to the bolt and torque according to the relevant instructions. When using a Tube, place friction pad between foot and Tube.

## ALIGNMENT INSTRUCTIONS

### ***Bench Alignment (Fig. 3)***

#### **Alignment Goal**

Alignment reference line (B) should:

- pass through midpoint of socket at either the patellar tendon level or at the ischial tuberosity level (D).
- fall at the 1/3 mark on the inside of the foot cover.

**Note:** Prioritize knee alignment over foot alignment if there is a mismatch.

#### **Alignment Instructions**

1. Position the foot so that the alignment reference line (B) falls at the 1/3 mark on the inside of the foot cover (with the foot cover and shoe on). Consider the external rotation of the foot.
2. Use the applicable adapters to connect either the socket or the knee to the foot and establish the correct knee center height.
3. If using a prosthetic knee: Position knee according to knee alignment instructions
4. On the lateral side of the socket, make a first mark (alignment reference point) at the midpoint of the socket at either patellar tendon level or the ischial tuberosity level (D). Make a second mark at the midpoint of the socket distally (E). Draw a line through both marks.

5. Position the socket so the alignment reference line (B) goes through the alignment reference point (D).
6. Set the correct socket angles for flexion/extension and abduction/adduction.
7. If using a prosthetic knee: Use the applicable adapters to connect the socket to the prosthesis.

### **Static Alignment**

- Make sure the patient stands with equal weight on both legs.
- Check for correct prosthesis length.
- Check internal / external rotation.
- Check for correct load on toe and heel.

### **Dynamic Alignment**

Ensure the patient is familiar with the functioning of the device.

The heel to toe action can be influenced by:

- Heel Stiffness
- Anterior-Posterior positioning of device
- Dorsi-Plantarflexion
- Shoe Characteristics

Consider the following actions if needed:

#### **Symptoms**

- Device comes to flat position too early (patient feels he/she is sinking into a hole)
- Climbing over the toe requires extra energy
- Toe feels too stiff
- Knee hyperextends

#### **Action**

- Check Wedge instructions (if applicable)
- Shift socket anterior (or device posterior)
- Consider dorsiflexion
- Check heel of the shoe and shoe performance

#### **Symptoms**

- Rapid heel to toe movement
- Poor control over prosthesis at initial contact
- Minimal energy return feeling
- Too little push off from the toe
- Knee becomes unstable

#### **Action**

- Check Wedge instructions (if applicable)
- Shift socket posterior (or device anterior)
- Consider plantarflexion
- Check heel of the shoe and shoe performance

### **Heel Wedges (Fig. 4)**

The small, medium and large Wedges are used to change heel stiffness.

The Wedges can be trimmed using sharp scissors to customize stiffness.

For temporary Wedge placement, use tape to secure the wedge in position.

For permanent Wedge placement

- Roughen the upper and lower surface of the heel with abrasive paper.
- Apply instant adhesive on the lower side of the Wedge only.
- Locate in the foot/heel junction and position before adhesive sets.
- If the foot has a split toe, remove a thin slice in the middle by cutting with a sharp knife through the split in the carbon blades.

For removal, the adhesive may be softened by soaking in acetone or cyanoacrylate adhesive remover.

## **USAGE**

### ***Cleaning and care***

Clean with a damp cloth and a mild soap. Dry with a cloth after cleaning.

### ***Environmental Conditions***

The device is Waterproof if used with Tube.

The device is Weatherproof if used with Pyramid.

A Waterproof device can be used in a wet or humid environment and submerged in up to 3-meter-deep water for a maximum of 1 hour.

It can tolerate contact with: Salt water, chlorinated water, perspiration, urine, and mild soaps.

It can also tolerate occasional exposure to sand, dust, and dirt.

Continuous exposure is not allowed.

A Weatherproof device can be used in a wet or humid environment and can tolerate being splashed by fresh water (e.g., rain), no submersion is allowed.

No contact with salt water or chlorinated water is allowed.

Dry with a cloth after contact with fresh water or humidity.

Clean with fresh water after exposure to other liquids, chemicals, sand, dust, or dirt and dry with a cloth.

## **MAINTENANCE**

The device and the overall prosthesis should be examined by a healthcare professional. Interval should be determined based on patient activity.

### ***Noise from Foot***

Noise may occur if sand or debris is present in device. In that case, the healthcare professional should doff the foot, clean it with the help of compressed air and replace the Spectra sock if it is damaged.

## **REPORT OF SERIOUS INCIDENT**

Any serious incident in relation to the device must be reported to the manufacturer and relevant authorities.

## **DISPOSAL**

The device and packaging must be disposed of in accordance with respective local or national environmental regulations.

## **LIABILITY**

Össur prosthetic devices are designed and verified to be safe and compatible in combination with each other and custom-made prosthetic sockets with Össur adapters, and when used in accordance with their intended use.

Össur does not assume liability for the following:

- Device not maintained as instructed by the instructions for use.
- Device assembled with components from other manufacturers.
- Device used outside of recommended use condition, application, or environment.

## **Compliance**

This device has been tested according to ISO 10328 standard to two million load cycles.

Depending on patient activity this may correspond to 2-3 years of use.

ISO 10328 Label		
Category	Weight (Kg) Moderate Impact	Label Text
1	52	ISO 10328-P3*-52 kg
2	59	ISO 10328-P3*-59 kg
3	68	ISO 10328-P3-68 kg
4	77	ISO 10328-P3-77 kg
5	88	ISO 10328-P4-88 kg
6	100	ISO 10328-P5-100 kg
7	116	ISO 10328-P6-116 kg
8	130	ISO 10328-P7-130 kg
9	147	ISO 10328-P8-147 kg

\*Actual test loads reflect maximum body mass

ISO 10328 - "P" - "m"kg \*)

\*) Body mass limit not to be exceeded!

For specific conditions and limitations of use  
see manufacturer's written instructions on  
intended use!

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