

ASPIRE™ M3



Aspire M3 knee is a single-axis weight activated stance control knee. The adjustable swing phase allows for a more harmonious gait.

- Adjustable swing phase
- Adjustable brake sensitivity
- Adjustable extension assist



Patient Profile:	
Amputation Level:	Transfemoral
Impact Level:	Low to Moderate
Maximum Patient Weight:	100kg (220lbs)

Product Information:	
Build Height:	120mm
Weight of Knee:	0,445kg
Flexion:	150°



SINGLE AXIS



WEIGHT
ACTIVATED
LOCK



MECHANICAL
CONTROL




SINGLE SPEED



EXTENSION
ASSIST

CLEARANCE

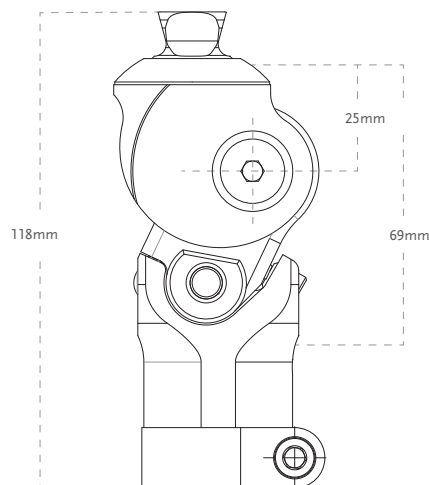

KITS

Part#	Description
1710513	Aspire M3