



MCPDriver®

A robust, custom, functional solution.

The MCPDriver is a body-driven prosthesis designed for finger and partial-hand difference at the proximal phalanx of digits two through five.

Adding an MCPDriver restores digit length to help recreate natural grip patterns. The MCPDriver excels at restoring precision and power grasps, including: pinch, key, and cylindrical grasps as well as grip stability.

The durable metal frame and robust components allow users to return to demanding lifestyles with confidence.

IDEAL ANATOMICAL PRESENTATIONS FOR FIT:

This device is driven by an intact MCP joint and suspension is maintained with enough residuum to engage the ring.

Our intake form will walk you through the steps of ordering a device using a Naked Prosthetics sizing kit, casting a positive hand mold, submitting detailed pictures of both hands, and a video showing range of motion.





OTHER COMFORT ASPECTS INCLUDE:

- Silicone interface of backplate adds comfort and cushion.
- Interchangeable rings allow for volume fluctuation.
- Adjustable washers can be moved to get possible alignment.
- Freedom of natural abduction and adduction allow the wearer to use the device.
- Conductive tip option for touchscreen use.

COLOURS + TEXTURES:

Location of texture may vary depending on order. Visit npdevices.com to see a variety of combinations.

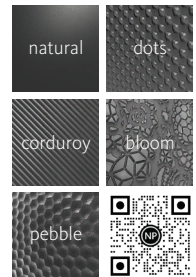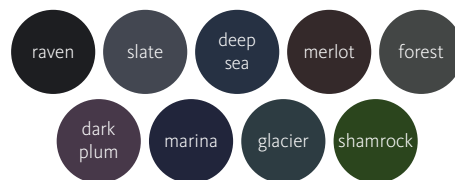

**Please note that due to the manufacturing process, variations in the colour, shade, and tone of nylon components and anodised fasteners may occur.*