



THE ALL NEW WATERPROOF

Proprio Foot®

Give your patients the best of both worlds – Bionic ankle motion and Pro-Flex carbon dynamics combine to normalise the walking experience as the world changes around them. Proprio Foot is waterproof and built for life. Our proven design has been updated with improved battery life and durability enhancements inside and out. It is fully submersible and carries IP67 ingress protection. All of this is now covered by a no hassle 36 month warranty.



The All New Waterproof Proprio Foot


Coupled with recent speed and adaptation enhancements, Proprio Foot is now suitable for a wide range of patient profiles; from the active hiker to the community ambulator.

- Motorised 33° ROM to pre-position foot. Foot blade provides 19.5° ROM in stance¹
- Adaptive ankle function reduces the energy cost of walking³
- Maintain full length toe lever, symmetry, and optimised carbon dynamics across a variety of surfaces⁴⁻⁷
- Active dorsiflexion to reduce the likelihood of trips and falls^{2,8}
- You work very hard optimising your patients' alignment, Auto-alignment helps them maintain it in the real world



- 1 IP67* water rating
- 2 New covers offer a sleek look and more durable finish
- 3 Updated twin screw design provides improved long term durability
- 4 Easy water drainage from prosthesis
- 5 Improved battery life – now up to 2 full days without charging
- 6 Streamlined setup and enhanced auto-adjustment with Össur Logic
- 7 Unity Vacuum Pump - also now waterproof

Product Information	
Amputation Level:	Transfemoral and Transtibial
Impact Level:	Low to Moderate
Maximum Patient Weight:	125kg
Waterproof Rating:	IP67*
Warranty:	3 Years standard warranty. Extendable up to 6 years of total coverage

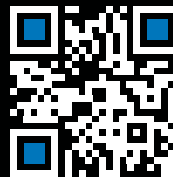

Scan or click here to learn more about the New Waterproof Proprio Foot



SANDAL TOE



HEEL HEIGHT ADJUSTABILITY



SPLIT TOE



ANATOMICAL FOOT



WATERPROOF*

*Proprio Foot will tolerate exposure to freshwater, salt water, perspiration, urine, mild soaps, sand, dust, and dirt. Exposure to chlorinated water is not recommended as it may compromise the actuator grease/lubrication, affecting functionality of the foot. Structural integrity is not affected in any way. Devices experiencing issues related to chlorine exposure can be serviced under warranty.

¹⁻⁸References can be viewed on the Proprio Foot page linked above.

FOLLOW ÖSSUR ON



WWW.OSSUR.CO.NZ

Össur New Zealand

TEL 0800 369 524

FAX 0800 448 265

customercare.nz@ossur.com

