



Pro-Flex Terra



Live Through Every Step

Hello,

Introducing the groundbreaking Pro-Flex Terra foot! This versatile foot, packed with innovation, is suitable for a wide range of users. From those with lower activity levels to high performers, Pro-Flex Terra fulfills the needs of your users as they change during the day.

A single foot is now sufficient for working in the office, going shopping, or running with friends, delivering performance broad enough to meet all the diverse daily requirements. Find out and discover all the ins- and outs of Pro-Flex Terra!

Pro-Flex Terra Features

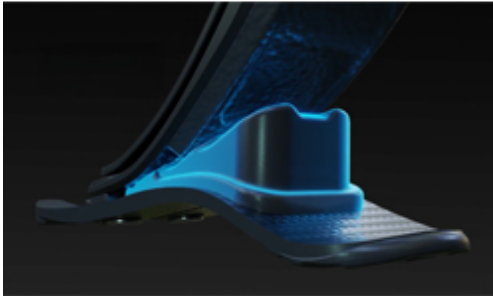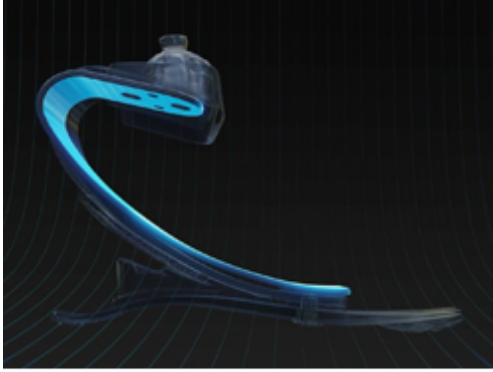

✓ 1. Pre-Tensioned Carbon

- Increased flexibility provides greater range of motion
- Increased comfort while standing
- Same force to achieve higher energy return

✓ 2. Single Anchor Point with Gliding Contact Pad

- Accommodates minor shifts in weight while standing
- Provides a neutral range during midstance by allowing the point of contact to move as weight shifts from heel to toe
- Provides multi-axial rotation for better ground adaptation

✓ 3. Pre-Compressed Foam

- Provides progressive stiffness across a range of impact levels
- Smoother transition from heel to toe
- Increased dampening at higher walking speeds

This is Pro-Flex Terra

Join us in revolutionising the industry together!



Upcoming Academy Courses:
Pro-Flex Family - Introducing the New Pro-Flex Terra

Dates: Wednesday 12 February, 2025

Wednesday 9 April, 2025

Location: Össur Academy facility in Melbourne

Join us for an in-depth, comprehensive session where you will gain valuable insights into the Pro-Flex family of feet, including their functions, applications, and alignment. Our expert-led training will empower you with the knowledge to enhance your understanding of these advanced products and their medical benefits. Upon completion, you will be able to:

- Differentiate between the various foot designs.
- Understand the medical benefits of each Pro-Flex foot, ensuring appropriate prescription for optimal patient outcomes.
- Recognise which feet are waterproof, weatherproof, and compatible with Unity systems.
- Accurately prescribe the correct foot design.

Register Today!

Resources

Please contact your clinical specialist for more information or to arrange a demonstration.

[Brochure](#)

[Catalogue Page](#)

[Mechanical Foot Comparison Chart](#)



Össur ANZ

Sydney: 26 Ross St, Parramatta, NSW 2151 | 1300 123 268

Melbourne: 9 Market Dr, Bayswater, VIC 3153 | 1300 123 268

Auckland: 6/80 Westpoint Dr, Hobsonville 0619 | 0800 369 524

www.ossur.com.au

© 2025 Össur

Want to change the types of emails you receive from Össur?

[Update your preferences?](#) | [Privacy Policy](#)